



Domus 2s

Alternating Pressure Redistribution System

- Suitable for patients with low to medium risk of pressure ulcers
- One piece mattress with slope design to relieve pressure on the heel
- Medical grade PVC pad for better hygiene and safety
- An ideal transition mattress for homecare




Domus 2s

Alternating Pressure Redistribution System



Anti-Sliding Extension Flap
Can be tucked under the standard mattress to prevent pad from sliding

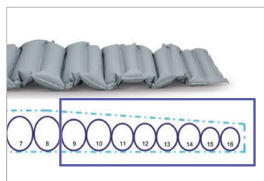


Heel Slope
Allow the pad to fill in the space behind the Achilles tendon to distribute pressure over the entire lower legs, reducing the risk of heel breakdown



Adjustable Hangers
Can secure control unit to fit most of beds



Weight Indicator Reference
To facilitate optimized pressure setting

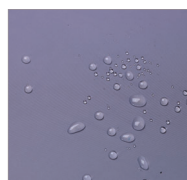

Nylon/PU Cover (Optional)
Vapor permeable and water resistant with basic flame retardant level



Water Resistance




Flame Resistance



Biocompatibility



Vapour Permeable

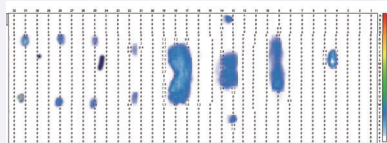

Pressure Relief During Alternating

Pressure mapping reveals that bladder A & B take turns to inflate and deflate can significantly relieve the body interface pressure

Specification

Pump	Dimension	9.8" x 4.9" x 3.3" / 25 x 11.2 x 9.5 cm
	Weight	2.65 lb / 1.2 Kg
	Case Material	Flame retardant ABS
	Supply Voltage	AC100-120V 60Hz, 0.1A
Mattress	Operating Cycle	10 min
	5" Overlay mattress	
	Dimension	78" x 33" x 5" : 5" for 1-8 cells, 4.5" for 9-10 cells, 4" for 11-14 cells, 3.5" for 15-16 cells
	Cell Height	16 x 5" ~ 3.5" cells
	Weight	6.83 lb / 3.1 Kg
	Material	Medical - grade PVC
	Maximum Patient Weight	260 lb / 120 Kg
	Flame Retardant Standards	CAL 117
	HCPCS code	E0372

Specifications are subject to change for improvement without notice

Ordering Information (For North America Only)

	Domus 2s
System	PP02101
Mattress	PM05301
Pump	PP02001